

VATT Spectrograph

STA0520A CCD
2688×512 pixels, 15 microns
ITL Serial Number 8228
Detector operating temperature -110C
Gain = 1.3 e/DN with setting 10
Gain = 2.6 e/DN with setting 5
For both gains the full well is ADC limited (65k DN), so this is 85,000 or 170,000 electrons.
Noise is 3.4 electrons (3.2 in overscan region)


Documentation

- [General Instructions](#)
- [Instructions on Guiding](#)

Maintenance

- Argon Calibration Lamp was replaced on 4/26/2013.
- Note: There is some concern the lamps are burning out faster than they should.

QE Curve

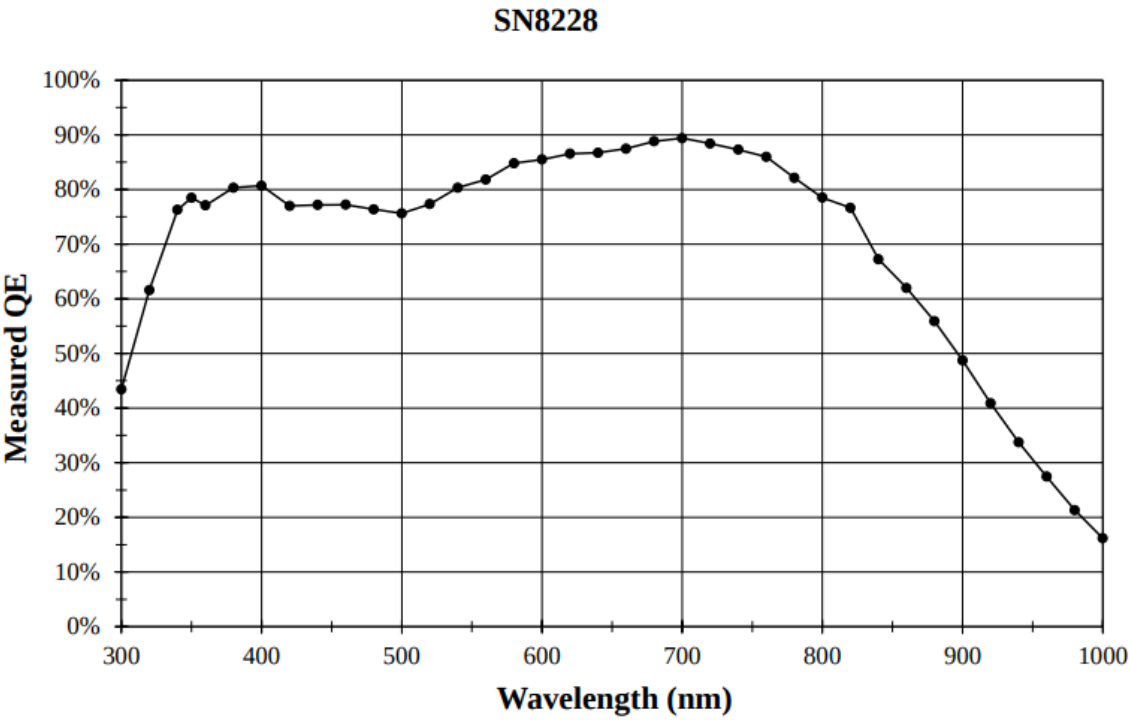

From:
<https://lavinia.as.arizona.edu/~tscopewiki/> - **MOON**

Permanent link:
<https://lavinia.as.arizona.edu/~tscopewiki/doku.php?id=vatt:vattspec&rev=1617248271>

Last update: **2021/03/31 20:37**

